

# **SIPOne**<sup>™</sup>

Suncombe Ltd, Jade House, Lockfield Avenue, Brimsdown, Enfield, Middlesex, EN37JY, United Kingdom *T*+44(0)20-8443-3454 *F*+44(0)20-8443-3969 *E* info@suncombe.com *W* www.suncombe.com





# About Us

Suncombe Ltd supply Cleaning and hygienic processing technology for the BioPharma, Medical, Food and other critical hygiene industries. From CIP to PartsWashers, we are uniquely capable of offering you a complete washing and processing solution.

Formed in 1961, we have a tremendous amount of experience in-house and provide a high level of technical and engineering expertise.

## Suncombe SIPOne<sup>™</sup> Range Details

The Suncombe SIPOne<sup>™</sup> has been developed as a stand-alone module for Sterilising In Place. Available in mobile and static configurations, the range has been developed over 45 years for providing a robust, reliable means of supplying repeatable, validateable, controllable SIP sterilisation. Stainless enclosures, with safety interlocked doors, to locate the equipment, can be supplied to contain the equipment for SIP, to ensure Health and Safety. Every unit is specifically designed for its particular application and is developed by Suncombe engineers together with the client personnel to provide the optimum sterilisation system and a simple validation trail. The range is based on standard modules with individual units custom designed and built - inhouse - for your specific application. Basic automation or 21CFR11 compliant reporting.

### **Key Features**

- 316L Stainless steel wetted parts, smooth, crevice-free construction
- PTFE elastomers and tri-clamp connections
- Enclosed head orbital welding to EN287/15614
- BioPharma instrumentation & Diaphragm valves
- Certified surface finishes to 0.1 ra and optionally electropolished




Typical Static SIPOne<sup>TM 46</sup>





Vessel in Enclosure for SIP



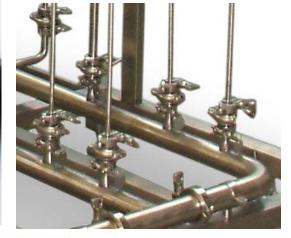




SIMATIC PANEL

### Automation

ControlSuite#2 and 3<sup>™</sup> - Automatic control system with f0 algorithm, temperature and time configurable parameters, I/O faceplates, Mimic diagrams, Hard wired interfaces, Printer, Repeatable automatic cycles, User control with passwords, In built troubleshooting and diagnostic ability, profibus, ethernet and hard wired interafces, Printer, 21CFR11 Datalogger and batch reports on ControlSuite#3.

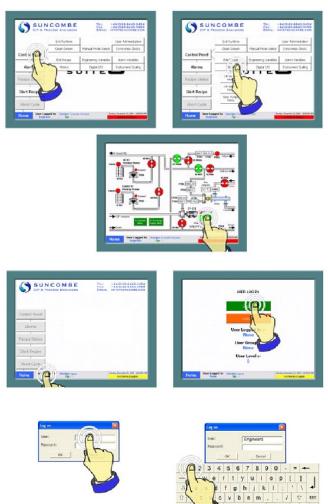







## **Support**

The Suncombe Technical Support and Customer Care departments' obligation is to provide total customer support. This support starts at the proposal stage, continues throughout the contract and thereafter. A personal customer care agent will be assigned to you at quotation stage and will support you throughout. The department currently supports all of our clients over the last 50 years. Support includes advise and surveys, installation and commissioning, User training, call-out repairs and service contracts for preventative maintenance.

